

MYCOTOXIN URINE COLLECTION INSTRUCTIONS

What's in the box?

- Test requisition form
- Urine transfer cup
- Vacuum Urine tube
- Green cap for shipping
- Absorbent Pad
- Biohazard Bag
- Ice Pak
- Prelabeled UPS Clinical Shipping Pak

Use the chart below for reference.
Urine's color should be between 2-5



**Discontinue binders
for 72 hours prior to
collection**



- Please allow 7-10 business days for results. Since this is a test ordered by your practitioner, you should request all results from their office.
- Urine must be medium to dark yellow, if you notice your urine is light colored or clear, please dispose, rinse out the collection cup, let it air out and recollect the next day.

TIP: discontinue liquids 2 hours before bed and avoid emptying your bladder in the middle of the night.

www.kbmodiagnostics.com



*Read both sides of these
instructions thoroughly
before starting the
collection process*

SCAN HERE!




For more detailed
collection and
shipping instructions

- 1 Collect the first morning urine into the urine cup, THEN transfer to the transport tube (**a,b,c**)
- 2 After transferring the sample into the tube, remove the white cap on the transport tube and replace it with the green cap before freezing and shipping
- 3 Fill out the patient's name, date of birth (DOB), and date of collection (DOC) on the transport tube label (c) fine tip permanent ink suggested
- 4 Freeze sample with sample tube upright and freeze ice pak laying flat, until completely frozen
- 5 Complete ALL highlighted parts of the test requisition form
- 6 Fold test requisition form and place in the outer pouch of the biohazard bag to protect it from ice pak condensation
- 7 Once ice pak and sample are completely frozen, your sample is ready for shipping
- 8 Seal the ice pack, sample tube, and absorbent pad in the biohazard bag's Ziploc section. Place everything with the paperwork in the original kit box
- 9 Place the box in the UPS clinical pak and make sure is completely sealed
- 10 Go to www.ups.com to search for the nearest drop-off spot. Drop your sample at any authorized UPS drop off location

a



b

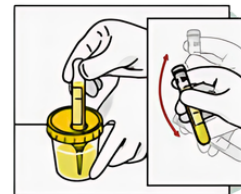

c

